

Composite Wetting Resin

Description and Indications:

Composite Wetting Resin is a 45% filled, radiopaque light cured wetting resin. It may be used during the placement of composite if the composite has become dry during incremental layering, or if the oxygen inhibition layer has been removed or disturbed. Composite Wetting Resin helps to improve the glide of the instrument during sculpting and contouring.

Procedure:

1. Remove luer lock caps and attach the Inspiral® Brush tip to syringe.
2. Apply a thin layer to cured composite increment as needed. Air or vacuum thin (vacuum excess resin using 18-20 gauge Endo-Eze® tip on Luer Vacuum adapter).
3. If using Composite Wetting Resin for instrument glide, apply a small quantity to instrument or surface of soft composite.

Precautions:

1. Any resin can be sensitizing to skin or can cause allergic reaction if individual is already sensitized. Therefore, wash area of skin exposed to resin completely to remove the material. If dermatitis, rash or any other allergic reaction is observed, see a physician.
2. When dispensing intraorally, never force syringe plungers. If there is excess resistance to flow, replace tip with a new one and verify flow before proceeding.
3. Use only recommended tips with each specific material. Tips are disposable.
4. We recommend using Ultradent® Syringe Covers sealed with the Ultradent® Impulse Sealer for a protective barrier to prevent cross-contamination.

Composite Wetting Resin



D-Composite Wetting Resin/F-Agent de mouillage pour composite/NL-Composietbevestigingshars/
IT-Resina per bagnare i compositi/E-Resina Hidratante de Composite/P-Resina Humificadora para
Compósito/S-Vättningsresin för komposit/DK-Komposit Befugtningresin/FIN-Yhdistelmämuovin
kostutusresiini/GR-Σύνθετη ρητίνη για διαβροχή



Inspir® Brush Tip



English

Apply a thin layer to cured composite increment as needed. Air or vacuum thin (vacuum excess resin using 18-20 gauge Endo-Eze® tip on Luer Vacuum adapter).

If using Composite Wetting Resin for instrument glide, apply a small quantity to instrument or surface of soft composite.


German

Eine dünne Schicht auf ausgehärtete Komposit-Schicht nach Bedarf aufbringen. Mit Luft oder durch Absaugen ausdünnen (Überschüssiges Resin mit Hilfe von 18-20 gauge Endo-Eze-Kanülen und Vacuum-Adapter absaugen).

Wenn Composite Wetting Resin zur Instrumentenbenetzung verwendet werden soll, geringe Menge auf Instrument oder auf die Oberfläche des weichen Komposits auftragen.

Ultradent syringes have an expiration date stamped on the side of the syringe consisting of one letter and three numbers. The letter is a lot number used for manufacturing purposes and the three numbers are the expiration date. The first two numbers are the month, and the third number is the last number of the year.

Composite Wetting Resin

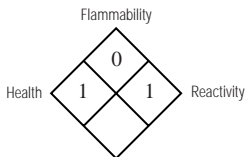

HAZARD RATING

- 4 = Severe
- 3 = Serious
- 2 = Moderate
- 1 = Slight
- 0 = Minimal

MSDS Included.

Shelf life: 24 months. Store at room temperature.

U.S. federal law restricts this device to sale by or on the order of a dentist.

For professional use only. Keep out of reach of children.

For immediate reorder and/or complete descriptions of Ultradent's product line, refer to Ultradent's catalog or call Toll Free **1-800-552-5512**. Outside U.S. call (801) 572-4200.

© Copyright 2003 Ultradent Products, Inc. All Rights Reserved.

U.S. Patent Nos.: 4,997,371; 5,269,684; 5,534,562; and 5,708,052. International Patents Granted and Pending.

ULTRADENT
PRODUCTS, INC.

Manufactured by: Ultradent Products, Inc.
505 West 10200 South, South Jordan, Utah 84095

#23887.2 100903